

User Manual MK-dent LED BU8012EM


Technical Data

Art.Nr. BU8012EM
Luminance 3,0V - 3,5V (best light with 3,2V)


Installation



The MK-dent LED works with direct current, that means the LED has two different poles. If the LED **does not work** after installation, you need to re/insert the LED after rotating it 180°. To achieve the optimal luminance, we recommend a voltage of 3,2 V.

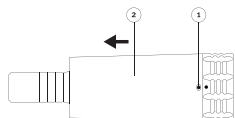

☞ Push button 1 and pull off sleeve 2 socket and take it out.

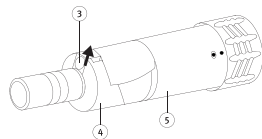

- ☞ Press the porter 4 to the body 5.
- ☞ Push the LED 3 with a screw driver in direction of arrow out of the socket. Take it out.
- ☞ Place the new LED in the socket. Pay attention to the correct position of contacts.
- ☞ Push the LED into the socket and reattach sleeve.
- ☞ Turn LED 180°, in case it doesn't light.



WARNING

WARNING

Risk of injury by high heat addition

Do not touch the LED during use.